

SOCOREX®

ACURA® ELECTRO PIPETTES



DWK LIFE SCIENCES
DURAN WHEATON KIMBLE

SOCOREX®

Acura® *electro* Pipettes



Your Electronic Choice

The Acura® *electro* pipettes are the latest upgrade of a successful line that makes electronic pipetting more versatile, simpler and safer than ever.

The Acura *electro* pipettes offer a unique cost-saving advantage. Once the pipette initial package is purchased, subsequent purchases consist of the appropriate volumetric module or pipette alone. The initial package contains the charging stand that holds up to three pipettes or battery packs.

The microprocessor controlled pipette features intuitive software, eliminating any complicated programming steps. There are six working modes that are easy to access and many volumes to select from. If you are looking for unsurpassed metrological performance and stability, choose the Acura *electro* pipettes.

Product Advantages:

- Optimal ergonomics, lightweight
- Easy, intuitive, self-teaching program
- Large display, reversible for right or left reading
- Swift-interchange battery pack
- Fast charge, long working autonomy
- Adjustable tip ejector* fitting most tips
- Pipetting cycle counter
- 27 interchangeable volumetric modules - all fitting same control unit
- Two-year warranty

*Socorex patented



Product Advantages and Instrument Performance



Micropipettes

- Reduced shaft length and conical end improve driveability
- Easy access to microtubes and microplates
- Optional longer volumetric modules
- Simple conversion by user at any time
- Unsurpassed metrological performance and stability

Macropipettes

- Interchangeable nozzle protection filter
- Adapters for Pasteur pipette and straw tips
- Outstanding accuracy and precision

Multichannel

- Large selection of 8 and 12-channel models
- Lightweight and perfect hand fitting
- Sequential tip ejection
- Also fits any single channel volumetric module

Ordering Information

Initial package includes pipette, charging stand and accessories. Additional pipettes (alone) supplied with individual QC certificate, pipette tips samples and operating instructions.

Cat. No.	Cat. No.	Volume	Division	Inaccuracy (E%)			Imprecision (CV%)			Tip Style
Initial Package	Pipette Alone			Min. vol.	Mid. Vol.	Max. Vol.	Min. Vol.	Mid. Vol.	Max. Vol.	
Micropipettes Acura® electro 926XS										
W870902-A	W870902	0.1 - 2µL	0.01µL	<± 2.5 % ¹⁾	<± 1.2 %	<± 0.9 %	< 2.5 % ¹⁾	< 1.5 %	< 0.8 %	Ultra 10µL
W870904-A	W870904	0.5 - 10µL	0.05µL	<± 1.2 % ²⁾	<± 0.8 %	<± 0.6 %	< 1.5 % ²⁾	< 0.7 %	< 0.35 %	Ultra 10µL
W870906-A	W870906	0.5 - 10µL	0.05µL	<± 1.2 % ²⁾	<± 0.8 %	<± 0.6 %	< 1.7 % ²⁾	< 0.8 %	< 0.4 %	200µL
W870908-A	W870908	1 - 20µL	0.1µL	<± 1.2 % ²⁾	<± 0.6 %	<± 0.5 %	< 1.2 % ²⁾	< 0.4 %	< 0.3 %	200µL
W870910-A	W870910	2.5 - 50µL	0.25µL	<± 1.0 % ²⁾	<± 0.6 %	<± 0.5 %	< 0.7 % ²⁾	< 0.3 %	< 0.25 %	200µL
W870912-A	W870912	5 - 100µL	0.5µL	<± 1.0 % ²⁾	<± 0.6 %	<± 0.5 %	< 0.7 % ²⁾	< 0.3 %	< 0.2 %	200µL
W870914-A	W870914	10 - 200µL	1.0µL	<± 1.0 % ²⁾	<± 0.6 %	<± 0.4 %	< 0.6 % ²⁾	< 0.2 %	< 0.15 %	200µL
W870916-A	W870916	50 - 1000µL	5.0µL	<± 0.8 % ²⁾	<± 0.5 %	<± 0.4 %	< 0.4 % ²⁾	< 0.15 %	< 0.1 %	1000µL
Macropipettes Acura® electro 936										
W870918-A	W870918	0.1 - 2mL	0.01mL	<± 1.5 % ²⁾	<± 1.0 %	<± 0.5 %	< 0.6 % ²⁾	< 0.3 %	< 0.15 %	2mL
W870920-A	W870920	0.25 - 5mL	0.05mL	<± 1.2 % ²⁾	<± 0.8 %	<± 0.5 %	< 0.6 % ²⁾	< 0.3 %	< 0.15 %	5mL
W870922-A	W870922	0.5 - 10mL	0.05mL	<± 1.0 % ²⁾	<± 0.7 %	<± 0.5 %	< 0.5 % ²⁾	< 0.2 %	< 0.15 %	10mL
8-Channel Pipettes Acura® electro 956										
W870924-A	W870924	0.5 - 10µL	0.05µL	<± 3.5 % ²⁾	<± 1.5 %	<± 1.0 %	< 3.0 % ²⁾	< 0.9 %	< 0.7 %	Ultra 10µL
W870926-A	W870926	2.5 - 50µL	0.25µL	<± 1.0 % ²⁾	<± 0.9 %	<± 0.8 %	< 1.0 % ²⁾	< 0.6 %	< 0.4 %	200µL
W870928-A	W870928	10 - 200µL	1.0µL	<± 0.9 % ²⁾	<± 0.7 %	<± 0.6 %	< 0.6 % ²⁾	< 0.4 %	< 0.25 %	200µL
W870930-A	W870930	20 - 350µL	5.0µL	<± 1.0 % ²⁾	<± 0.8 %	<± 0.6 %	< 0.6 % ²⁾	< 0.4 %	< 0.25 %	350µL
12-Channel Pipettes Acura® electro 956										
W870932-A	W870932	0.5 - 10µL	0.05µL	<± 3.5 % ²⁾	<± 1.5 %	<± 1.0 %	< 3.0 % ²⁾	< 0.9 %	< 0.7 %	Ultra 10µL
W870934-A	W870934	2.5 - 50µL	0.25µL	<± 1.0 % ²⁾	<± 0.9 %	<± 0.8 %	< 1.0 % ²⁾	< 0.6 %	< 0.4 %	200µL
W870936-A	W870936	10 - 200µL	1.0µL	<± 0.9 % ²⁾	<± 0.7 %	<± 0.6 %	< 0.6 % ²⁾	< 0.4 %	< 0.25 %	200µL
W870938-A	W870938	20 - 350µL	5.0µL	<± 1.0 % ²⁾	<± 0.8 %	<± 0.6 %	< 0.6 % ²⁾	< 0.4 %	< 0.25 %	350µL

Performance values obtained in Forward mode with double-distilled water at constant temperature (± 0.5°C) comprised between 20 and 25°C in accordance with ISO 8655.

Performance measured at: ¹⁾ 0.5 µL, ²⁾ 10% of nominal volume

SOCOREX® Acura® *electro* Electronic Pipette

Great Built-in Flexibility

Scroll through mode selection by pressing a single button.



Forward Mode

Aspiration and pipetting of set volume. Suitable for all applications.



Reverse Mode

Aspiration in excess, followed by dosing of set volume. Enhances reproducibility below 20µL. Especially recommended for viscous and foaming liquids.



Stepper Mode

Tip filling and step-by-step distribution. Best suited to aliquot samples.



Dilution Mode

Aspiration of 2 or 3 different volumes for restitution in one shot. An easy way to dilute samples.



Tactile Mode

"Start and stop" liquid measurement, titration and gel loading by simple touch of activation button.



Mixing

"Up and down" liquid flow in the tip / vial.



1

Natural Ergonomics 1

The Acura® *electro* offers shape, balance and working ergonomics resembling those of a manual instrument. But comfort is greatly improved, thus pipetting and result consistency are outstanding.



2

Unique "Right or Left" Display Reading 2

Display converts instantly to right or left side reading. All information is clearly presented and visible at a glance during programming, pipetting and calibration steps.



3

Variable Working Speed 3


Located on the front side, speed selector allows swift change even during the pipetting process. In addition, slow speed access at any time by gentle touch of start button.



4

Pipetting Cycle Counter 4

A simple double click gives access to the number of pipetting cycles performed since last zeroing.

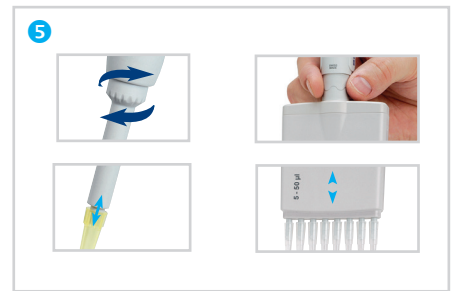


5

Tip Ejection - Easier Than Ever 5

Maximal efficiency obtained with a large surface, ergonomically located ejector button. The shaft height adjustment system called Justip™ (4mm span) - controlled by efficient click-stops - allows a wide selection of tips to tightly fit the nozzle.

The multichannel bowed shape ejector head permits effortless, sequential tip ejection.

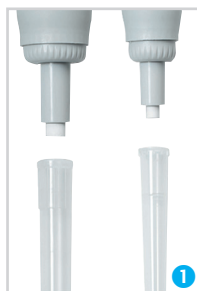

1

Nozzle Protection Filter 1


Macro models accommodate a filter protecting instrument against liquid overflow and contamination. Located in the nozzle, it can be removed and changed easily.



2

Pasteur Pipette - Straw Tip Adapters 2

Optional accessories fit 2mL glass Pasteur pipettes or 5mL straw tip instead of Polypropylene tips.

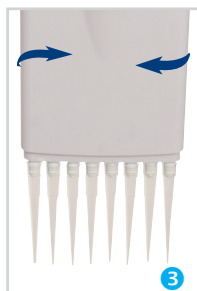

3

Optimal Working Position 3

360° rotation of multichannel volumetric modules allows selection of the best working position.



4

Time Efficient Maintenance 4

Pipette construction limits maintenance to a minimum - no tool required to remove volumetric module. The pipetting cycle counter facilitates maintenance monitoring. Whenever sterilization is needed, volumetric modules can be autoclaved at 121°C/250°F.



5

Initial Package 5

Recommended when purchasing an Acura® *electro* for the first time, each initial package contains one electronic pipette, QC certificate and operating instructions, charging stand, power supply, additional battery pack, pipette tips samples. Thereafter, additional pipettes can be purchased alone and charged on existing stand.

SOCOREX® Acura® *electro* Electronic Pipette

Interchangeable Volumetric Modules



Calibration 1

Control software gives immediate access to calibration menu. Performance check possible over two or three independent volumes (V_{min} , V_{mid} and V_{max}). New settings entered directly from the keyboard. Error message warns of any inaccurate plunger movement.

Compact Charging Rack 2

Space saving stand holds up to three battery packs simultaneously. The ideal charging and storage accessory for spare batteries, available at any time.


Pipette Charging Stand 3

Stand allows charging and storing up to three electronic pipettes or battery packs simultaneously. Each position is fitted with a red / green charging LED. One stand supplied with every initial package.

Long Working Autonomy 4

NiMH battery pack exchangeable instantly. Fast battery charging (<1.5 hours) and extended working autonomy (>3000 consecutive pipetting cycles). Battery charge level clearly visible on display. Automatic energy saving stand-by when unused.

One control unit fits a choice of 27 volumetric assemblies available separately. This greatly extends working possibilities and makes electronic pipetting affordable to all budgets. Disassembling and re-assembling without any tool contributes towards maximal flexibility. Calibration specificities of extra modules, as set by factory QC or by pipette owner, are retained in instrument memory.

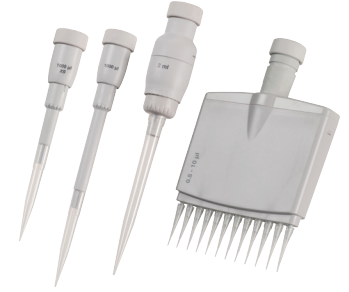

Volumetric Modules

Cat. No.	Volume Range	Tip Style
Single Channel Microvolumes Reduced Shaft Length		
W870012	0.1 - 2 μ L	Ultra 10 μ L
W870014	0.5 - 10 μ L	Ultra 10 μ L
W870016	0.5 - 10 μ L	200 μ L
W870018	1 - 20 μ L	200 μ L
W870020	2.5 - 50 μ L	200 μ L
W870022	5 - 100 μ L	200 μ L
W870024	10 - 200 μ L	200 μ L
W870026	50 - 1000 μ L	1000 μ L
Regular Shaft Length		
W870112	0.1 - 2 μ L	Ultra 10 μ L
W870114	0.5 - 10 μ L	Ultra 10 μ L
W870116	0.5 - 10 μ L	200 μ L
W870118	1 - 20 μ L	200 μ L
W870120	2.5 - 50 μ L	200 μ L
W870122	5 - 100 μ L	200 μ L
W870124	10 - 200 μ L	200 μ L
W870126	50 - 1000 μ L	1000 μ L
Single Channel Macrovolumes		
W870128	0.1 - 2mL	2mL
W870130	0.25 - 5mL	5mL
W870132	0.5 - 10mL	10mL
8-Channel Microvolumes		
W870134	0.5 - 10 μ L	Ultra 10 μ L
W870136	2.5 - 50 μ L	200 μ L
W870138	10 - 200 μ L	200 μ L
W870140	20 - 350 μ L	350 μ L
12-Channel Microvolumes		
W870142	0.5 - 10 μ L	Ultra 10 μ L
W870144	2.5 - 50 μ L	200 μ L
W870146	10 - 200 μ L	200 μ L
W870148	20 - 350 μ L	350 μ L

Accessories

Cat. No.	Description	Qty/Pack
W851342	Nozzle Protection Filters for 2 and 5mL Models, PP Fibers	250
851346	Nozzle protection Filters for 10mL Model, PP Fibers	100
W820023	Pasteur Pipette Adapter Nozzle for 5mL Model	1
W820124	Straw Tip Adapter for 5mL Models	1

Charging Accessories

Cat. No.	Description	Qty/Pack
W820019	Pipette Charging Stand, 3 Charging Positions, Holds Pipettes or Spare Batteries	1
W820020	Compact Charging Stand, 3 Charging Positions, Holds Spare Batteries Only	1
W870001	Spare Battery Pack (NiMH, 4.8V)	1
W870002	Spare Battery Duo Pack (NiMH, 4.8V)	2
W820022-A	Power Supply, 100-240V	1

SOCOREX® Acura® *electro* Electronic Pipette

Interchangeable Volumetric Modules



Volume Range
926 XS micro
0.1 - 2µL
0.5 - 10µL
0.5 - 10µL
1 - 20µL
2.5 - 50µL
5 - 100µL
10 - 200µL
50 - 1000µL
936 macro
0.1 - 2mL
0.25 - 5mL
0.5 - 10mL
956 multi
0.5 - 10µL
2.5 - 50µL
10 - 200µL
20 - 350µL







Excellence in your hands

WBRO_062 03/2019



Authorized Distributor
www.netascientific.com
(800) 343-6015
orders@netascientific.com

DWK Life Sciences
1501 North 10th Street
Millville, NJ 08332
United States

US & Canada: 800.225.1437
Int'l: 856.825.1100
Fax: 856.825.1368
www.dwk.com